


Magnetic Bumper Set

ME-9885A

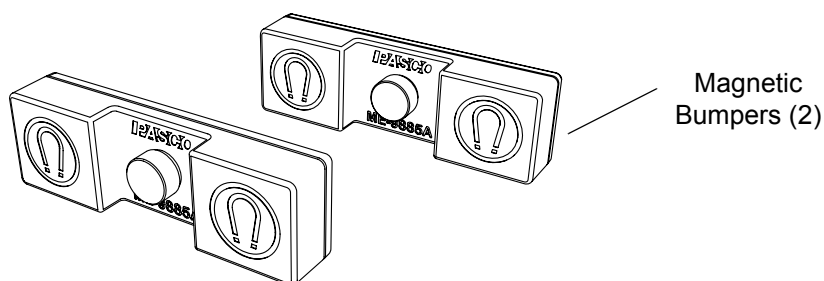

Included Equipment	Part Number
Magnetic Bumper (2)	003-09567

Additional Equipment Recommended	
Force Accessory Bracket	CI-6545
Force Sensors	CI-6537, CI-6746 or PS-2104
PASCO carts	ME-6951, ME-6950 or ME-9454

Introduction

The Magnetic Bumper Set includes a pair of strong magnetic bumpers. These bumpers are ideal for measuring forces during elastic collisions between carts. Simply remove the hooks on the Force Sensors and replace with the magnetic bumpers. The bumper has a metal post which assures that it is mounted in the correct orientation.

Suggested Activities

Elastic Collisions

1. Remove the hooks from both Force Sensors
2. Screw a Magnetic Bumper into each of the sensors, Note the location of the alignment post on the bumpers and the hole on the sensor face.
3. Connect the Force Sensors to PASCO carts or Force Accessory Bracket



Limited Warranty

For a description of the product warranty, see the PASCO catalog.

Contacting Technical Support

Address: PASCO scientific
10101 Foothills Blvd.
Roseville, CA 95747-7100

Web: www.pasco.com

Tech Support: support@pasco.com

Phone: (916) 786-3800
(800)772-8700

Copyright Notice

The PASCO scientific 012-09743A Magnetic Bumper Set Instruction Manual is copyrighted with all rights reserved.